

The Vesta-Syde SQ Ready to Use (RTU) products are a family of phosphate-free, one-step cleaner and disinfectants formulated for use on hard, non-porous surfaces. These Ready to Use formulations are based on the Vesta-Syde SQ Disinfectant Concentrate, and are available in dilutions of 1.6% v/v (1:64) and 0.8% v/v (1:128) to kill a broad spectrum of common bacteria and fungi. These two dilutions are available in both non-sterile and sterile, double-bagged packaging. The Vesta-Syde SQ RTU Disinfectants have an extensive documentation package to meet your validation objectives and are supported by the STERIS Technical Services team for individualized application and validation assistance.



## Vesta-Syde™ SQ Ready to Use Disinfectants

Easy to use, global quaternary ammonium disinfectants

### Features

Formulation is ready to use without any mixing

Broad spectrum efficacy against bacteria, fungi, and viruses using AOAC and EN protocols

Same key ingredients as Vesta-Syde SQ Disinfectant concentrate

Low residue formulation that is easily rinsed if required

Superior material compatibility

Active ingredients notified according to the Biocide Product Directive and registered with the US EPA

Sterile version is gamma irradiated and each lot USP sterility tested

Manufactured in accordance with cGMP with complete lot traceability

### Benefits

Easy to use and saves the time and labor of disinfectant preparation  
Eliminates errors in dilution preparation and the need to validate use procedures

Confirms efficacy against many common microbes

Simplifies validation requirements and offers multiple options for preparation

Ability to maintain visibly clean facilities with minimal rinsing

Allows for high frequency use on most surfaces


Available in most countries worldwide to facilitate program harmonization

Allows easy entry into controlled environments

Supports quality assurance and quality control goals

### Applications

- Surface disinfection and microbial control of:
  - » Floors
  - » Walls
  - » Ceilings
  - » Countertops
  - » BSL hoods
  - » Carts and other equipment
- Products may be applied by:
  - » Mopping
  - » Wiping
  - » Spraying

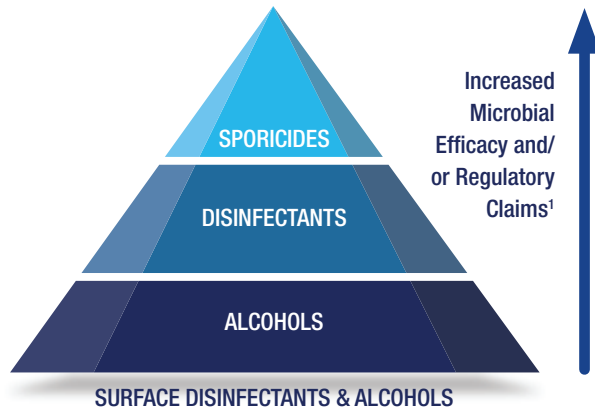

<sup>1</sup> Products that fall into the categories at the bottom of the pyramid are most frequently used and are generally not sporicidal. Progression up the pyramid indicates stronger performance overall and a broader spectrum of claims.

The foundation of any contamination control program is the disinfectant that should be used daily. Disinfectants should be used on a scheduled rotation with a sporicide based on your environmental data. Historically, multiple disinfectants have been rotated, but if your disinfectant has broad spectrum efficacy this may not be necessary. Vesta-Syde SQ RTU Disinfectants are easy to use while allowing for the use of Vesta-Syde SQ Disinfectant concentrate for flexibility in your contamination control program. Our Technical Service group can help customize your program to take advantage of these excellent products.

## Ordering Information (Worldwide)

### Vesta-Syde SQ 64 Ready to Use Disinfectant

Filtered for consistent quality for use in a variety of areas

Description		Europe	Canada
22 oz (651 mL) bottle (12 per case)	634277	634277EC	634277CA
32 oz (946 mL) bottle (12 per case)	634296	634296EC	634296CA
1 gallon (3.78 L) bottle (4 per case)	634208	634208EC	634208CA

### Vesta-Syde SQ 64 st Ready to Use Disinfectant

Double bagged and gamma irradiated through a validated process for easy introduction into your aseptic areas

Description		Europe	Canada
16 oz (473 mL) bottle (12 per case)	1S24Y4	1S24Y4EC	1S24Y4CA
23 oz (680 mL) bottle (12 per case)	1S24Y2	1S24Y2EC	1S24Y2CA
1 gallon bottle (3.78 L) (4 per case)	1S2308	1S2308EC	1S2308CA

### Vesta-Syde SQ 128 Ready to Use Disinfectant

Filtered for consistent quality for use in a variety of areas

Description		Europe	Canada
22 oz (651 mL) bottle (12 per case)	634377	634377EC	634377CA
32 oz (946 mL) bottle (12 per case)	634396	634396EC	634396CA
1 gallon bottle (3.78 L) (4 per case)	634308	634308EC	634308CA

### Vesta-Syde SQ 128 st Ready to Use Disinfectant

Double bagged and gamma irradiated through a validated process for easy introduction into your aseptic area

Description		Europe	Canada
16 oz (473 mL) bottle (12 per case)	1S25Y4	1S25Y4EC	1S25Y4CA
23 oz (680 mL) bottle (12 per case)	1S25Y2	1S25Y2EC	1S25Y2CA
1 gallon bottle (3.78 L) (4 per case)	1S2508	1S2508EC	1S2508CA